

# Anabolic Lab™

Analytical Lab Testing for  
Anabolic Steroid Harm Reduction

<https://anaboliclab.com>

Lab Report – November 18, 2016

Pharmacom Labs  
PHARMA Nan P100 (100 mg/mL)  
Nandrolone Phenylpropionate



If you like the attached report, please visit [our website](https://anaboliclab.com) to find out how you can support independent and objective lab testing.



Simec AG

Gärtnererweg 4 B CH-4665 Oftringen  
Telefon +41 62 752 83 08 Fax +41 62 752 83 09  
E-Mail: info@simec.ch www.simec.ch

# Test report

Number **C 48771**

**Object of analysis:** Pharmacom NAN P100

**Customer:** AnabolicLab.com

**Samples:** 1

**Batch:** 31j0c2m5vg

**Receipt:** 29.09.2016

**Subject:** Content, Identification Nandrolone Phenylpropionate

**Method(s):** HPLC, SOP M 006

| Sample / Analysis:                         | Method    | Status | Result      | Specification | complies |
|--|-----------|--------|-------------|---------------|----------|
| <b><u>Nandrolone Phenylpropionate:</u></b> |           |        |             |               |          |
| Identification                             | HPLC      | N      | complies    | complies      | yes      |
| Content                                    | HPLC      | N      | 103.1 mg/mL | 100 mg/mL     | -        |
| TAMC                                       | SOP M 006 | N      | < 1 CFU/mL  | -             | -        |
| TYMC                                       | SOP M 006 | N      | < 1 CFU/mL  | -             | -        |

Method-status: G = GMP A = accredited V = generally validated P = validated on product N = not validated  
E = external lab

**Remarks:** The sample(s) mentioned above have been analysed as they have been sent to us by the client. We have not controlled external processes like manufacturing, labelling, sampling, shipping, storage. These results are for information only and do not compensate for correct quality control by the manufacturer/distributor. Generally, pharmaceutical products have to be produced and distributed under full GMP/GDP regime, including GMP-compliant analyses with validated methods. It is not allowed to use this analysis in the context of doping/sports/competitions or any other illegal action, neither by the athlete nor by the tutor/trainer or any other person.

**Signatures:** controlled & released:

**Staff:**

**completed:**

18.11.2016


The enclosed results refer exclusively to the object of examination described above. The measuring accuracy is available on request. Only handwritten signatures are valid. Without written permission, it is not allowed to publish single parts of this report.

page 1 of 1

Report steroid analysis

HPLC: Nr. 6  
 Analysendatum: 17.11.2016 09:25:13  
 Vial: 60  
 Probenidentifikation: C48771 A  
 Operator: admin  
 Chrom.-Beding. ACN / Wasser, gradient  
 Säule Nr. 135, 100-5 C18 250 x 4.6 mm  
 Detektion: UV 240nm

Druckdatum: 18.11.2016 08:11:45  
 Injektvolumen: 20

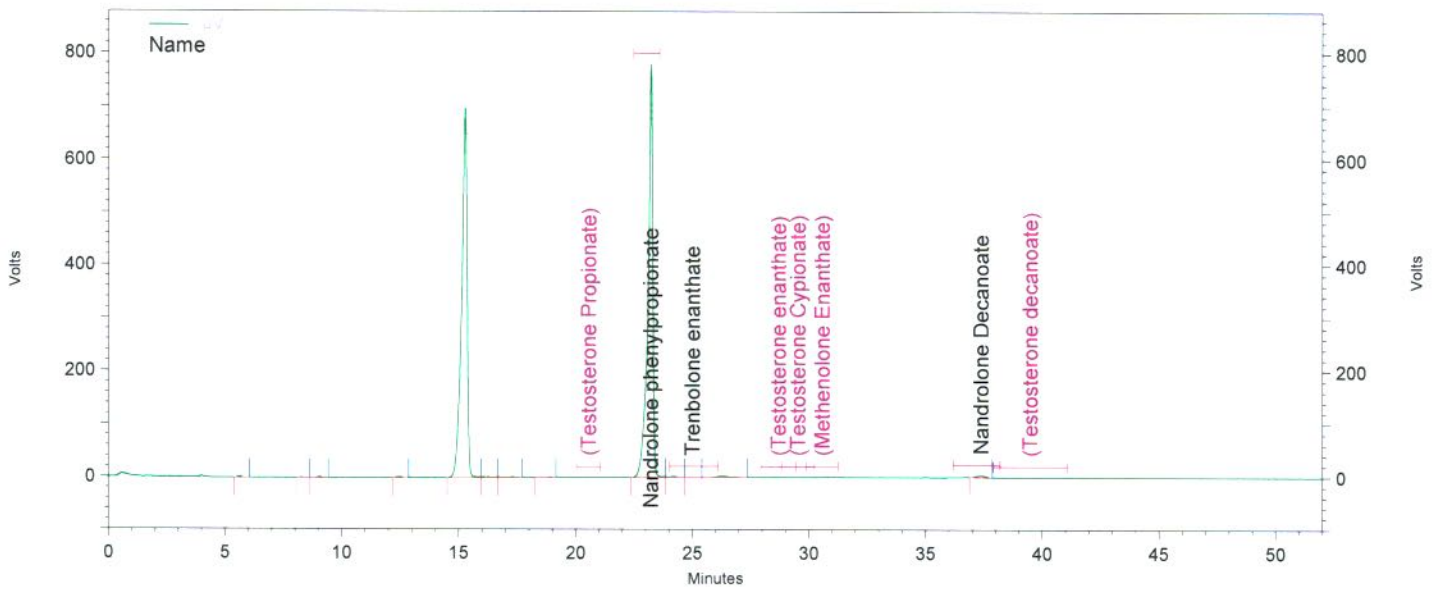

Channel A

UV Results

| Name                        | Time   | Area     | Amount    | Area Percent |
|-----------------------------|--------|----------|-----------|--------------|
|                             | 5.613  | 70646    | 0.000     | 0.077        |
|                             | 8.277  | 40510    | 0.000     | 0.044        |
|                             | 9.030  | 74725    | 0.000     | 0.081        |
|                             | 12.473 | 81044    | 0.000     | 0.088        |
|                             | 15.267 | 45337046 | 0.000     | 49.338       |
|                             | 16.263 | 115667   | 0.000     | 0.126        |
|                             | 17.303 | 100863   | 0.000     | 0.110        |
|                             | 18.913 | 47379    | 0.000     | 0.052        |
| Testosterone Propionate     |        |          | 0.000 BDL |              |
| Nandrolone phenylpropionate | 23.230 | 45231694 | 0.000     | 49.223       |
| Trenbolone enanthate        | 24.230 | 152703   | 0.000     | 0.166        |
|                             | 24.997 | 40543    | 0.000     | 0.044        |
| Testosterone enanthate      | 26.280 | 281044   | 0.000     | 0.306        |
| Testosterone Cypionate      |        |          | 0.000 BDL |              |
|                             |        |          | 0.000 BDL |              |


|              |        |        |       |       |
|--------------|--------|--------|-------|-------|
| Methenolone  |        |        | 0.000 | BDL   |
| Enanthate    |        |        |       |       |
| Nandrolone   | 37.377 | 317724 | 0.000 | 0.346 |
| Decanoate    |        |        |       |       |
| Testosterone |        |        | 0.000 | BDL   |
| decanoate    |        |        |       |       |

|        |  |          |       |         |
|--------|--|----------|-------|---------|
| Totals |  | 91891588 | 0.000 | 100.000 |
|--------|--|----------|-------|---------|

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\November 2016\09-11-2016\Diverse  
Steroide\020\_C48771 A-Rep1.dat  
Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Steroide A.met

Report steroid analysis

HPLC: Nr. 6  
 Analysendatum: 17.11.2016 10:18:32 Druckdatum: 18.11.2016 08:11:47  
 Vial: 60 Injektvolumen: 20  
 Probenidentifikation: C48771 A  
 Operator: admin  
 Chrom.-Beding. ACN / Wasser, gradient  
 Säule Nr. 135, 100-5 C18 250 x 4.6 mm  
 Detektion: UV 240nm

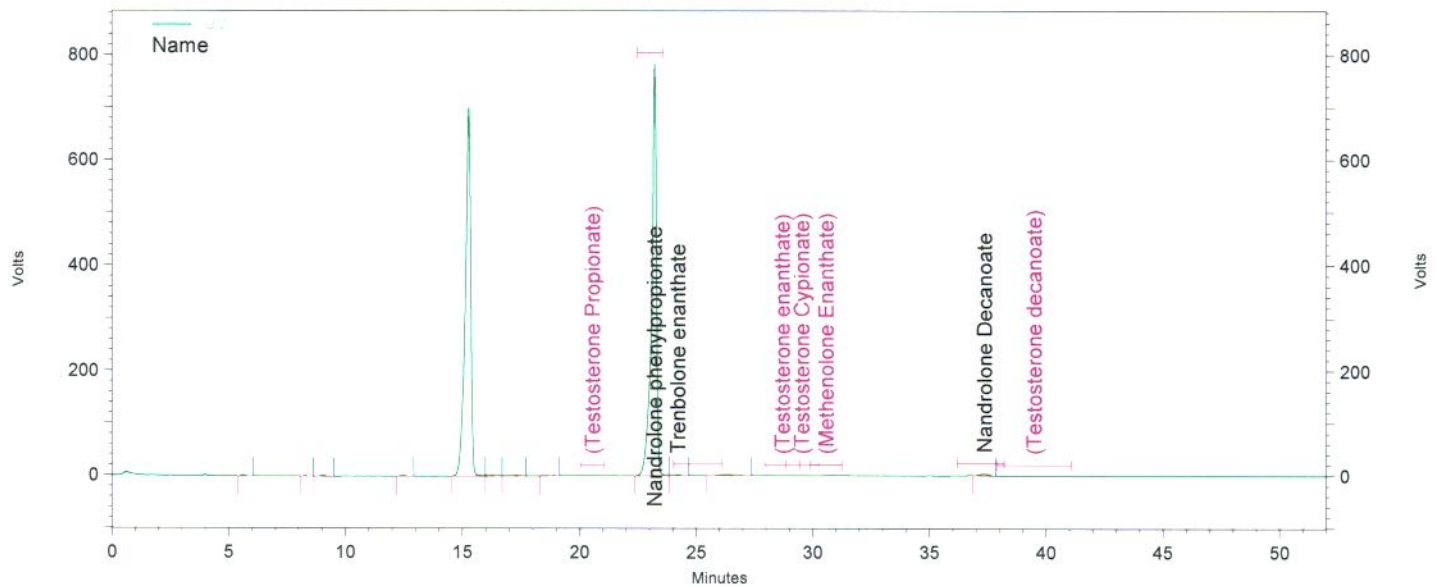

Channel A

| UV Results                  |        |          |           |              |
|-----------------------------|--------|----------|-----------|--------------|
| Name                        | Time   | Area     | Amount    | Area Percent |
|                             | 5.613  | 71509    | 0.000     | 0.077        |
|                             | 8.277  | 40243    | 0.000     | 0.044        |
|                             | 9.033  | 75362    | 0.000     | 0.082        |
|                             | 12.473 | 82167    | 0.000     | 0.089        |
|                             | 15.267 | 45562323 | 0.000     | 49.316       |
|                             | 16.263 | 119049   | 0.000     | 0.129        |
|                             | 17.307 | 99634    | 0.000     | 0.108        |
|                             | 18.910 | 47543    | 0.000     | 0.051        |
| Testosterone Propionate     |        |          | 0.000 BDL |              |
| Nandrolone phenylpropionate | 23.227 | 45528598 | 0.000     | 49.279       |
| Trenbolone enanthate        | 24.227 | 157761   | 0.000     | 0.171        |
|                             | 26.277 | 282084   | 0.000     | 0.305        |
| Testosterone enanthate      |        |          | 0.000 BDL |              |
| Testosterone Cypionate      |        |          | 0.000 BDL |              |
| Methenolone Enanthate       |        |          | 0.000 BDL |              |



|                           |        |        |           |       |
|---------------------------|--------|--------|-----------|-------|
| Nandrolone<br>Decanoate   | 37.360 | 322649 | 0.000     | 0.349 |
| Testosterone<br>decanoate |        |        | 0.000 BDL |       |


|        |  |          |       |         |
|--------|--|----------|-------|---------|
| Totals |  | 92388922 | 0.000 | 100.000 |
|--------|--|----------|-------|---------|

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\November 2016\09-11-2016\Diverse  
Steroide\020\_C48771 A-Rep2.dat  
Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Steroide A.met

Report steroid analysis

HPLC: Nr. 6  
 Analysendatum: 17.11.2016 11:11:52  
 Vial: 60  
 Probenidentifikation: C48771 A  
 Operator: admin  
 Chrom.-Beding.: ACN / Wasser, gradient  
 Säule: Nr. 135, 100-5 C18 250 x 4.6 mm  
 Detektion: UV 240nm

Druckdatum: 18.11.2016 08:11:50  
 Injektvolumen: 20

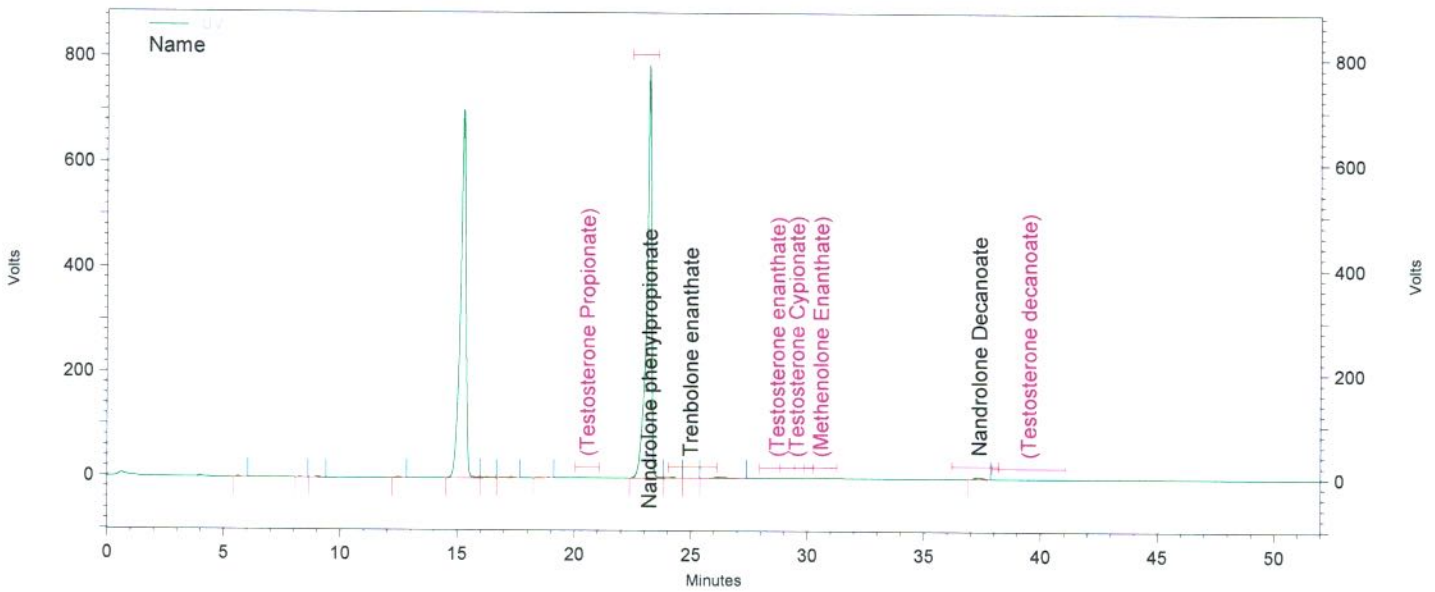

Channel A

UV Results

| Name                        | Time   | Area     | Amount    | Area Percent |
|-----------------------------|--------|----------|-----------|--------------|
|                             | 5.617  | 70974    | 0.000     | 0.077        |
|                             | 8.280  | 40471    | 0.000     | 0.044        |
|                             | 9.033  | 75066    | 0.000     | 0.081        |
|                             | 12.477 | 78514    | 0.000     | 0.085        |
|                             | 15.263 | 45636093 | 0.000     | 49.300       |
|                             | 16.263 | 114674   | 0.000     | 0.124        |
|                             | 17.300 | 98574    | 0.000     | 0.106        |
|                             | 18.907 | 46979    | 0.000     | 0.051        |
| Testosterone Propionate     |        |          | 0.000 BDL |              |
| Nandrolone phenylpropionate | 23.223 | 45606324 | 0.000     | 49.268       |
| Trenbolone enanthate        | 24.227 | 153757   | 0.000     | 0.166        |
|                             | 24.993 | 41055    | 0.000     | 0.044        |
| Testosterone enanthate      | 26.277 | 282350   | 0.000     | 0.305        |
| Testosterone Cypionate      |        |          | 0.000 BDL |              |

|              |        |        |           |       |
|--------------|--------|--------|-----------|-------|
| Methenolone  |        |        | 0.000 BDL |       |
| Enanthate    |        |        |           |       |
| Nandrolone   | 37.360 | 322511 | 0.000     | 0.348 |
| Decanoate    |        |        |           |       |
| Testosterone |        |        | 0.000 BDL |       |
| decanoate    |        |        |           |       |

|        |  |          |       |         |
|--------|--|----------|-------|---------|
| Totals |  | 92567342 | 0.000 | 100.000 |
|--------|--|----------|-------|---------|

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\November 2016\09-11-2016\Diverse  
Steroide\020\_C48771 A-Rep3.dat  
Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Steroide A.met






Report steroid analysis

HPLC: Nr. 6  
 Analysendatum: 17.11.2016 12:05:18  
 Vial: 61  
 Probenidentifikation: C48771 B  
 Operator: admin  
 Chrom.-Beding. ACN / Wasser, gradient  
 Säule Nr. 135, 100-5 C18 250 x 4.6 mm  
 Detektion: UV 240nm

Druckdatum: 18.11.2016 08:11:52  
 Injektvolumen: 20

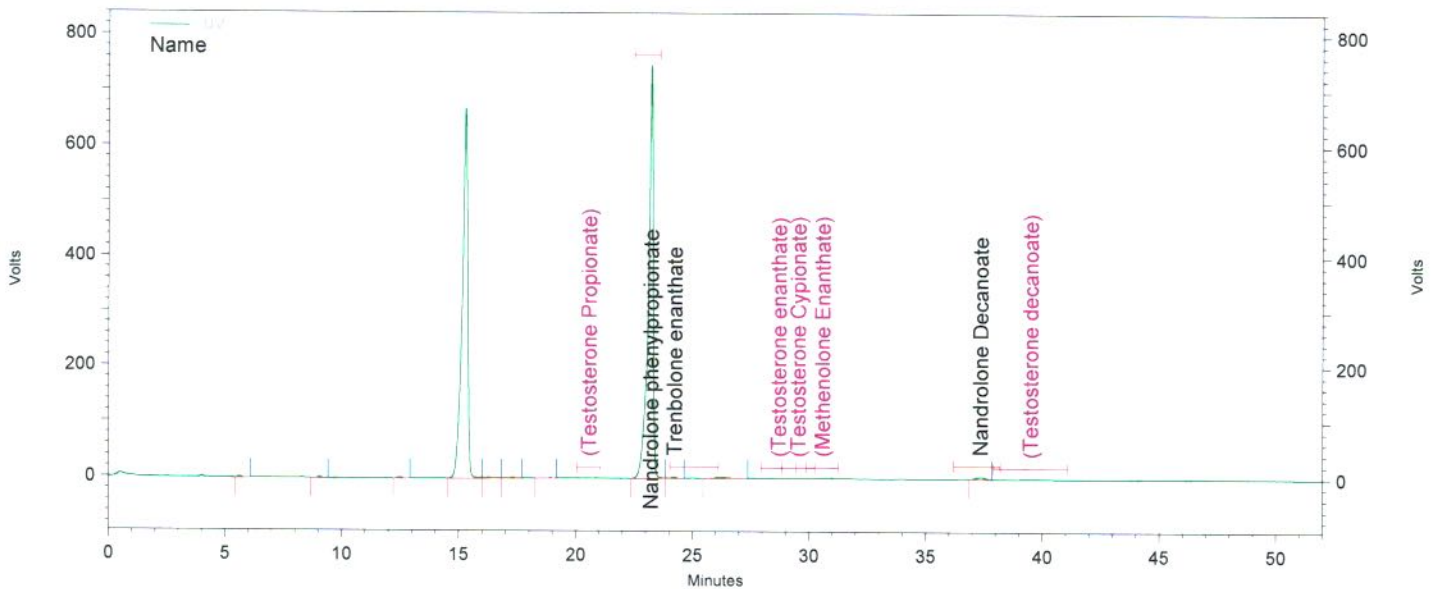

Channel A

UV Results

| Name                        | Time   | Area     | Amount    | Area Percent |
|-----------------------------|--------|----------|-----------|--------------|
|                             | 5.623  | 67255    | 0.000     | 0.076        |
|                             | 9.047  | 70785    | 0.000     | 0.080        |
|                             | 12.487 | 74766    | 0.000     | 0.085        |
|                             | 15.273 | 43383618 | 0.000     | 49.310       |
|                             | 16.273 | 118198   | 0.000     | 0.134        |
|                             | 17.307 | 90791    | 0.000     | 0.103        |
|                             | 18.913 | 47242    | 0.000     | 0.054        |
| Testosterone Propionate     |        |          | 0.000 BDL |              |
| Nandrolone phenylpropionate | 23.223 | 43368227 | 0.000     | 49.293       |
| Trenbolone enanthate        | 24.230 | 148943   | 0.000     | 0.169        |
|                             | 26.267 | 286210   | 0.000     | 0.325        |
| Testosterone enanthate      |        |          | 0.000 BDL |              |
| Testosterone Cypionate      |        |          | 0.000 BDL |              |
| Methenolone Enanthate       |        |          | 0.000 BDL |              |


|                           |        |        |           |       |
|---------------------------|--------|--------|-----------|-------|
| Nandrolone<br>Decanoate   | 37.350 | 324696 | 0.000     | 0.369 |
| Testosterone<br>decanoate |        |        | 0.000 BDL |       |

|        |  |          |       |         |
|--------|--|----------|-------|---------|
| Totals |  | 87980731 | 0.000 | 100.000 |
|--------|--|----------|-------|---------|

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\November 2016\09-11-2016\Diverse  
Steroide\021\_C48771 B-Rep1.dat  
Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Steroide A.met

Report steroid analysis

HPLC: Nr. 6  
 Analysendatum: 17.11.2016 12:58:43 Druckdatum: 18.11.2016 08:11:54  
 Vial: 61 Injektvolumen: 20  
 Probenidentifikation: C48771 B  
 Operator: admin  
 Chrom.-Beding. ACN / Wasser, gradient  
 Säule Nr. 135, 100-5 C18 250 x 4.6 mm  
 Detektion: UV 240nm

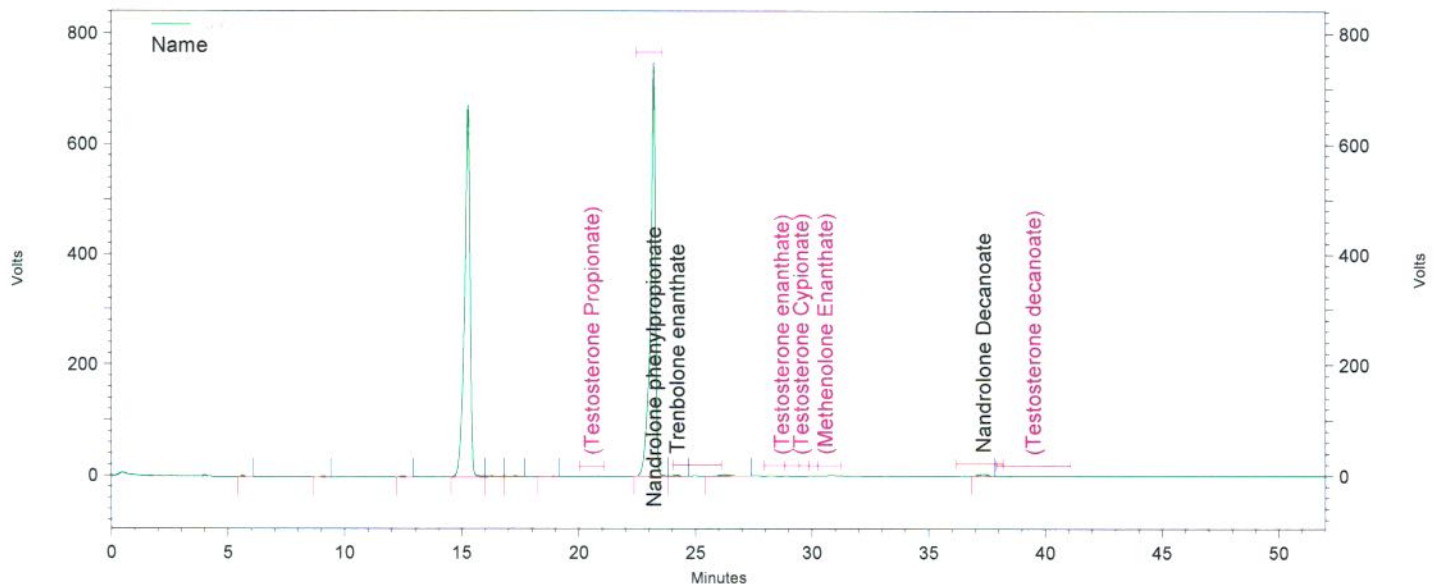

Channel A

| UV Results                  |        |          |           |              |
|-----------------------------|--------|----------|-----------|--------------|
| Name                        | Time   | Area     | Amount    | Area Percent |
|                             | 5.627  | 67634    | 0.000     | 0.077        |
|                             | 9.050  | 70680    | 0.000     | 0.080        |
|                             | 12.487 | 74616    | 0.000     | 0.085        |
|                             | 15.277 | 43492455 | 0.000     | 49.317       |
|                             | 16.277 | 113522   | 0.000     | 0.129        |
|                             | 17.310 | 91446    | 0.000     | 0.104        |
|                             | 18.913 | 45413    | 0.000     | 0.051        |
| Testosterone Propionate     |        |          | 0.000 BDL |              |
| Nandrolone phenylpropionate | 23.223 | 43461433 | 0.000     | 49.282       |
| Trenbolone enanthate        | 24.223 | 153996   | 0.000     | 0.175        |
|                             | 26.267 | 286274   | 0.000     | 0.325        |
| Testosterone enanthate      |        |          | 0.000 BDL |              |
| Testosterone Cypionate      |        |          | 0.000 BDL |              |
| Methenolone Enanthate       |        |          | 0.000 BDL |              |

|              |        |        |           |       |
|--------------|--------|--------|-----------|-------|
| Nandrolone   | 37.350 | 331907 | 0.000     | 0.376 |
| Decanoate    |        |        |           |       |
| Testosterone |        |        | 0.000 BDL |       |
| decanoate    |        |        |           |       |

|        |  |          |       |         |
|--------|--|----------|-------|---------|
| Totals |  | 88189376 | 0.000 | 100.000 |
|--------|--|----------|-------|---------|


Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\November 2016\09-11-2016\Diverse  
Steroide\021\_C48771 B-Rep2.dat  
Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Steroide A.met



Report steroid analysis

HPLC: Nr. 6  
 Analysendatum: 17.11.2016 13:52:07  
 Vial: 61  
 Probenidentifikation: C48771 B  
 Operator: admin  
 Chrom.-Beding. ACN / Wasser, gradient  
 Säule Nr. 135, 100-5 C18 250 x 4.6 mm  
 Detektion: UV 240nm

Druckdatum: 18.11.2016 08:11:56  
 Injektvolumen: 20

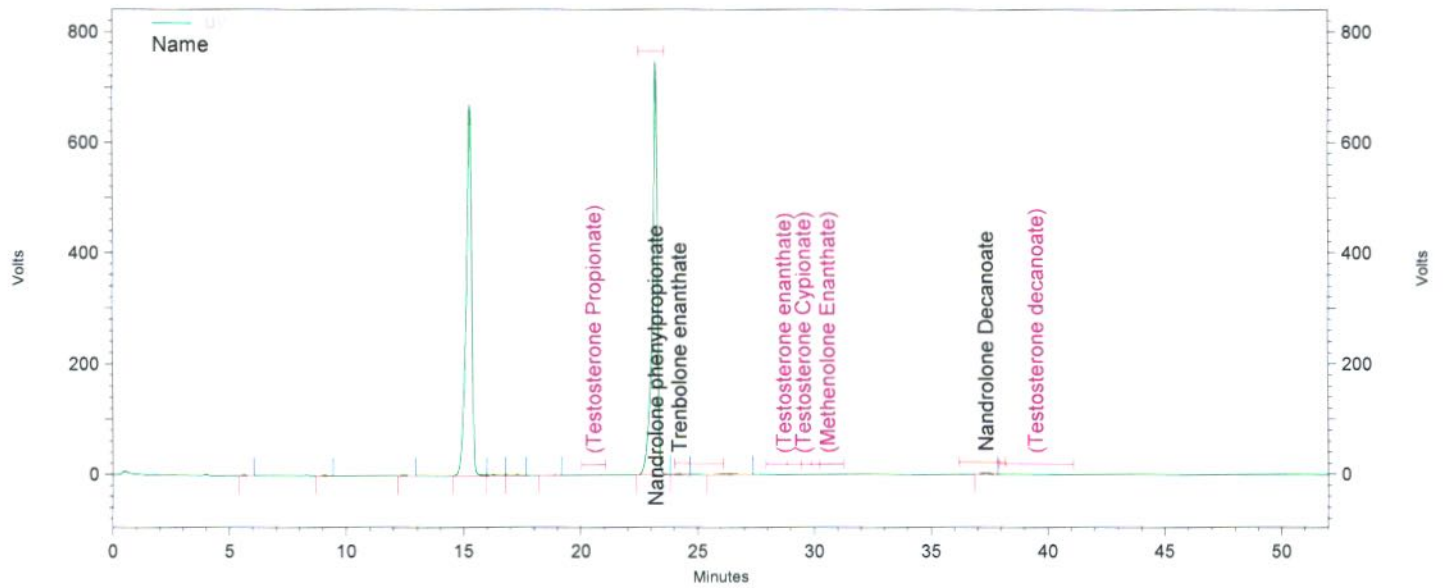

Channel A

UV Results

| Name                        | Time   | Area     | Amount    | Area Percent |
|-----------------------------|--------|----------|-----------|--------------|
|                             | 5.623  | 67246    | 0.000     | 0.076        |
|                             | 9.050  | 70887    | 0.000     | 0.080        |
|                             | 12.487 | 74398    | 0.000     | 0.084        |
|                             | 15.277 | 43522567 | 0.000     | 49.312       |
|                             | 16.273 | 115627   | 0.000     | 0.131        |
|                             | 17.313 | 90369    | 0.000     | 0.102        |
|                             | 18.910 | 46886    | 0.000     | 0.053        |
| Testosterone Propionate     |        |          | 0.000 BDL |              |
| Nandrolone phenylpropionate | 23.223 | 43509091 | 0.000     | 49.297       |
| Trenbolone enanthate        | 24.227 | 153902   | 0.000     | 0.174        |
|                             | 26.273 | 283256   | 0.000     | 0.321        |
| Testosterone enanthate      |        |          | 0.000 BDL |              |
| Testosterone Cypionate      |        |          | 0.000 BDL |              |
| Methenolone Enanthate       |        |          | 0.000 BDL |              |


|                           |        |        |           |       |
|---------------------------|--------|--------|-----------|-------|
| Nandrolone<br>Decanoate   | 37.353 | 324481 | 0.000     | 0.368 |
| Testosterone<br>decanoate |        |        | 0.000 BDL |       |

|        |  |          |       |         |
|--------|--|----------|-------|---------|
| Totals |  | 88258710 | 0.000 | 100.000 |
|--------|--|----------|-------|---------|

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\November 2016\09-11-2016\Diverse  
Steroide\021\_C48771 B-Rep3.dat  
Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Steroide A.met

Report steroid analysis

HPLC: Nr. 6  
 Analysendatum: 10.11.2016 00:36:35 Druckdatum: 14.11.2016 15:11:03  
 Vial: 3 Injektvolumen: 20  
 Probenidentifikation: STD Nandrolone Phenylprop und Test Enanthate  
 Operator: admin  
 Chrom.-Beding. ACN / Wasser, gradient  
 Säule Nr. 119, 100-5 C18 250 x 4.6 mm  
 Detektion: UV 240nm

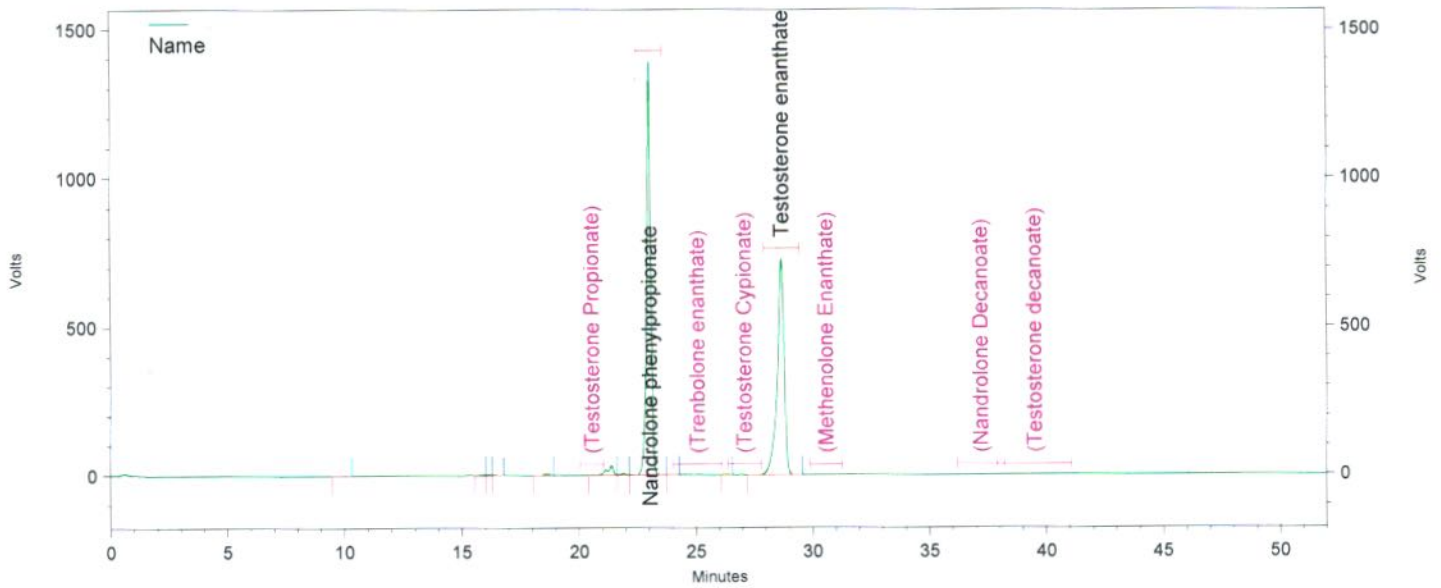

Channel A

UV Results

| Name                        | Time   | Area     | Amount    | Area Percent |
|-----------------------------|--------|----------|-----------|--------------|
|                             | 10.023 | 56231    | 0.000     | 0.043        |
|                             | 15.880 | 110105   | 0.000     | 0.084        |
|                             | 16.153 | 124989   | 0.000     | 0.095        |
|                             | 16.383 | 100505   | 0.000     | 0.076        |
|                             | 18.623 | 243321   | 0.000     | 0.185        |
| Testosterone Propionate     |        |          | 0.000 BDL |              |
|                             | 21.393 | 2200107  | 0.000     | 1.674        |
|                             | 21.917 | 248799   | 0.000     | 0.189        |
| Nandrolone phenylpropionate | 23.043 | 68790852 | 0.000     | 52.334       |
|                             | 24.017 | 48801    | 0.000     | 0.037        |
| Trenbolone enanthate        |        |          | 0.000 BDL |              |
|                             | 26.350 | 46422    | 0.000     | 0.035        |
| Testosterone Cypionate      |        |          | 0.000 BDL |              |
| Testosterone enanthate      | 28.700 | 59476590 | 0.000     | 45.248       |
| Methenolone Enanthate       |        |          | 0.000 BDL |              |


Nandrolone 0.000 BDL  
Decanoate  
Testosterone 0.000 BDL  
decanoate

|        |  |           |       |         |
|--------|--|-----------|-------|---------|
| Totals |  | 131446722 | 0.000 | 100.000 |
|--------|--|-----------|-------|---------|

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\November 2016\09-11-2016\Diverse  
Steroide\005\_STD Nandrolone Phenylprop und Test Enanthate-Rep1.dat  
Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Steroide A.met

Report steroid analysis

HPLC: Nr. 6  
 Analysendatum: 10.11.2016 01:30:01 Druckdatum: 14.11.2016 15:11:05  
 Vial: 3 Injektvolumen: 20  
 Probenidentifikation: STD Nandrolone Phenylprop und Test Enanthate  
 Operator: admin  
 Chrom.-Beding. ACN / Wasser, gradient  
 Säule Nr. 119, 100-5 C18 250 x 4.6 mm  
 Detektion: UV 240nm

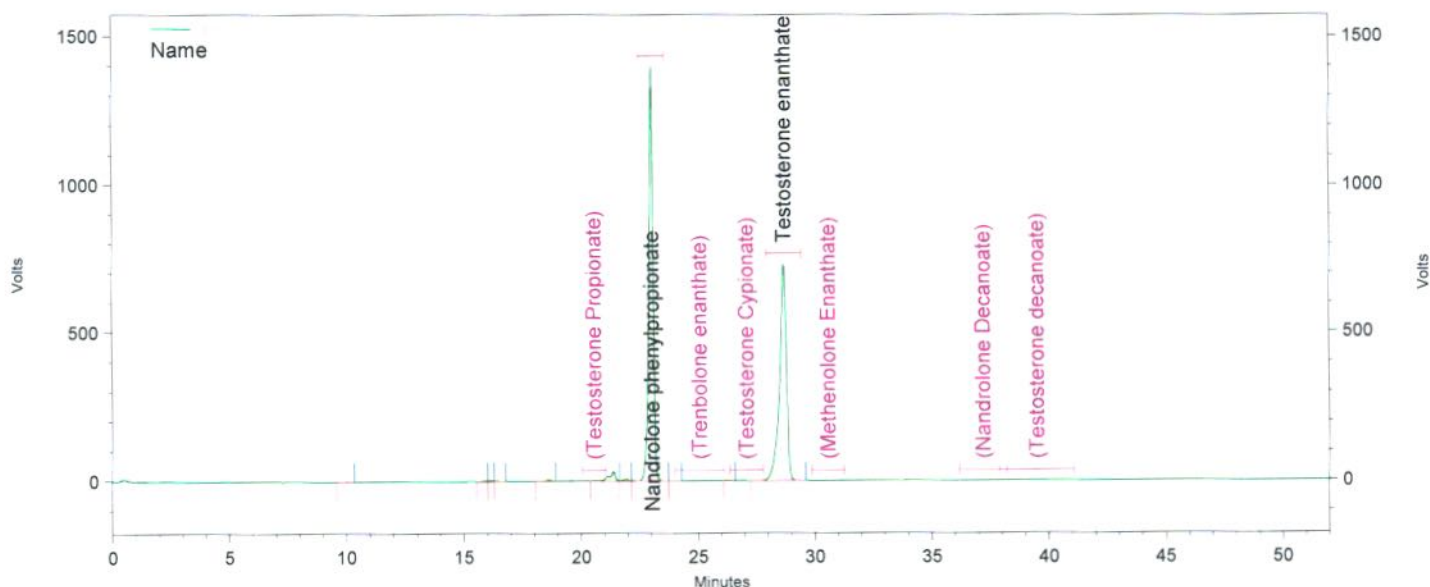

Channel A

| UV Results                  |        |          |           |              |
|-----------------------------|--------|----------|-----------|--------------|
| Name                        | Time   | Area     | Amount    | Area Percent |
|                             | 10.033 | 54836    | 0.000     | 0.042        |
|                             | 15.880 | 111350   | 0.000     | 0.085        |
|                             | 16.153 | 126021   | 0.000     | 0.096        |
|                             | 16.383 | 102212   | 0.000     | 0.078        |
|                             | 18.627 | 242430   | 0.000     | 0.184        |
| Testosterone Propionate     |        |          | 0.000 BDL |              |
|                             | 21.393 | 2208605  | 0.000     | 1.676        |
|                             | 21.917 | 251319   | 0.000     | 0.191        |
| Nandrolone phenylpropionate | 23.043 | 68934662 | 0.000     | 52.315       |
|                             | 24.020 | 48875    | 0.000     | 0.037        |
| Trenbolone enanthate        |        |          | 0.000 BDL |              |
|                             | 26.350 | 46747    | 0.000     | 0.035        |
| Testosterone Cypionate      |        |          | 0.000 BDL |              |
| Testosterone enanthate      | 28.700 | 59642563 | 0.000     | 45.263       |
| Methenolone Enanthate       |        |          | 0.000 BDL |              |




Nandrolone 0.000 BDL  
Decanoate  
Testosterone 0.000 BDL  
decanoate

|        |  |           |       |         |
|--------|--|-----------|-------|---------|
| Totals |  | 131769620 | 0.000 | 100.000 |
|--------|--|-----------|-------|---------|

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\November 2016\09-11-2016\Diverse  
Steroide\005\_STD Nandrolone Phenylprop und Test Enanthate-Rep2.dat  
Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Steroide A.met

Report steroid analysis

HPLC: Nr. 6  
 Analysendatum: 10.11.2016 02:23:23 Druckdatum: 14.11.2016 15:11:12  
 Vial: 3 Injektvolumen: 20  
 Probenidentifikation: STD Nandrolone Phenylprop und Test Enanthate  
 Operator: admin  
 Chrom.-Beding. ACN / Wasser, gradient  
 Säule Nr. 119, 100-5 C18 250 x 4.6 mm  
 Detektion: UV 240nm

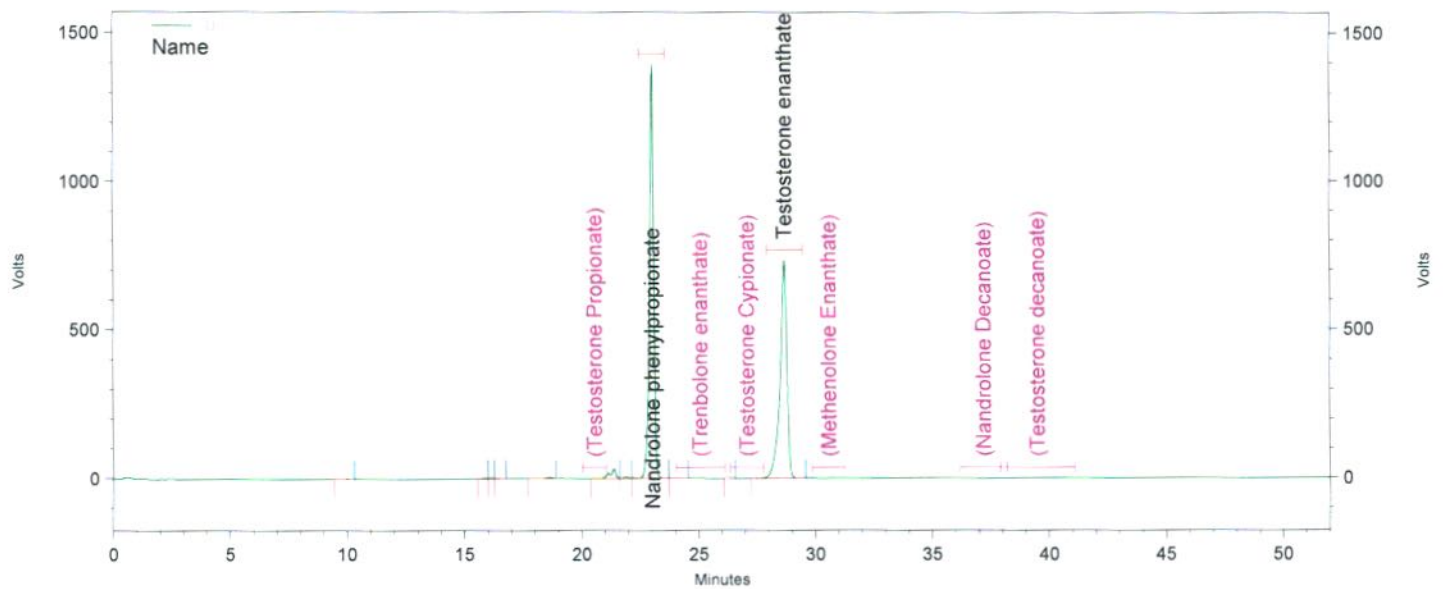

Channel A

UV Results

| Name                        | Time   | Area     | Amount    | Area Percent |
|-----------------------------|--------|----------|-----------|--------------|
|                             | 10.017 | 56354    | 0.000     | 0.043        |
|                             | 15.870 | 111927   | 0.000     | 0.085        |
|                             | 16.143 | 126701   | 0.000     | 0.096        |
|                             | 16.373 | 102349   | 0.000     | 0.077        |
|                             | 18.617 | 265971   | 0.000     | 0.201        |
| Testosterone Propionate     |        |          | 0.000 BDL |              |
|                             | 21.387 | 2218895  | 0.000     | 1.679        |
|                             | 21.907 | 260657   | 0.000     | 0.197        |
| Nandrolone phenylpropionate | 23.037 | 69127398 | 0.000     | 52.314       |
|                             | 24.013 | 81276    | 0.000     | 0.062        |
| Trenbolone enanthate        |        |          | 0.000 BDL |              |
|                             | 26.340 | 45973    | 0.000     | 0.035        |
| Testosterone Cypionate      |        |          | 0.000 BDL |              |
| Testosterone enanthate      | 28.690 | 59742433 | 0.000     | 45.211       |
| Methenolone Enanthate       |        |          | 0.000 BDL |              |


Nandrolone 0.000 BDL  
Decanoate  
Testosterone 0.000 BDL  
decanoate

|        |  |           |       |         |
|--------|--|-----------|-------|---------|
| Totals |  | 132139934 | 0.000 | 100.000 |
|--------|--|-----------|-------|---------|

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\November 2016\09-11-2016\Diverse  
Steroide\005\_STD Nandrolone Phenylprop und Test Enanthate-Rep3.dat  
Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Steroide A.met

Report steroid analysis

HPLC: Nr. 6  
 Analysendatum: 10.11.2016 03:16:46 Druckdatum: 14.11.2016 15:11:14  
 Vial: 3 Injektvolumen: 20  
 Probenidentifikation: STD Nandrolone Phenylprop und Test Enanthate  
 Operator: admin  
 Chrom.-Beding. ACN / Wasser, gradient  
 Säule Nr. 119, 100-5 C18 250 x 4.6 mm  
 Detektion: UV 240nm

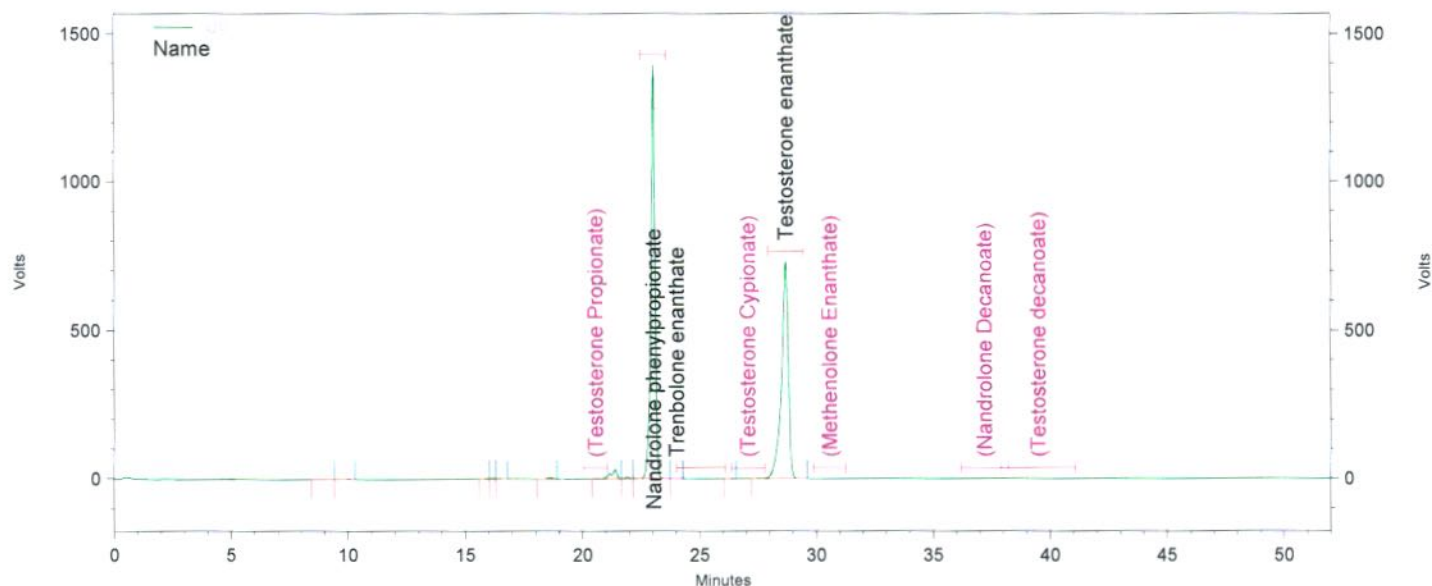

Channel A

UV Results

| Name                        | Time   | Area     | Amount    | Area Percent |
|-----------------------------|--------|----------|-----------|--------------|
|                             | 9.160  | 40643    | 0.000     | 0.031        |
|                             | 10.030 | 56145    | 0.000     | 0.043        |
|                             | 15.883 | 112157   | 0.000     | 0.085        |
|                             | 16.157 | 125518   | 0.000     | 0.095        |
|                             | 16.390 | 102263   | 0.000     | 0.077        |
|                             | 18.630 | 246083   | 0.000     | 0.186        |
| Testosterone Propionate     |        |          | 0.000 BDL |              |
|                             | 21.400 | 2211311  | 0.000     | 1.675        |
|                             | 21.923 | 250745   | 0.000     | 0.190        |
| Nandrolone phenylpropionate | 23.050 | 69060417 | 0.000     | 52.315       |
| Trenbolone enanthate        | 24.030 | 49641    | 0.000     | 0.038        |
|                             | 26.357 | 47484    | 0.000     | 0.036        |
| Testosterone Cypionate      |        |          | 0.000 BDL |              |
| Testosterone enanthate      | 28.707 | 59706968 | 0.000     | 45.229       |
| Methenolone Enanthate       |        |          | 0.000 BDL |              |

Nandrolone 0.000 BDL  
Decanoate  
Testosterone 0.000 BDL  
decanoate


|        |  |           |       |         |
|--------|--|-----------|-------|---------|
| Totals |  | 132009375 | 0.000 | 100.000 |
|--------|--|-----------|-------|---------|

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\November 2016\09-11-2016\Diverse  
Steroide\005\_STD Nandrolone Phenylprop und Test Enanthate-Rep4.dat  
Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Steroide A.met



Report steroid analysis

HPLC: Nr. 6  
 Analysendatum: 10.11.2016 04:10:07 Druckdatum: 14.11.2016 15:11:16  
 Vial: 3 Injektvolumen: 20  
 Probenidentifikation: STD Nandrolone Phenylprop und Test Enanthate  
 Operator: admin  
 Chrom.-Beding. ACN / Wasser, gradient  
 Säule Nr. 119, 100-5 C18 250 x 4.6 mm  
 Detektion: UV 240nm

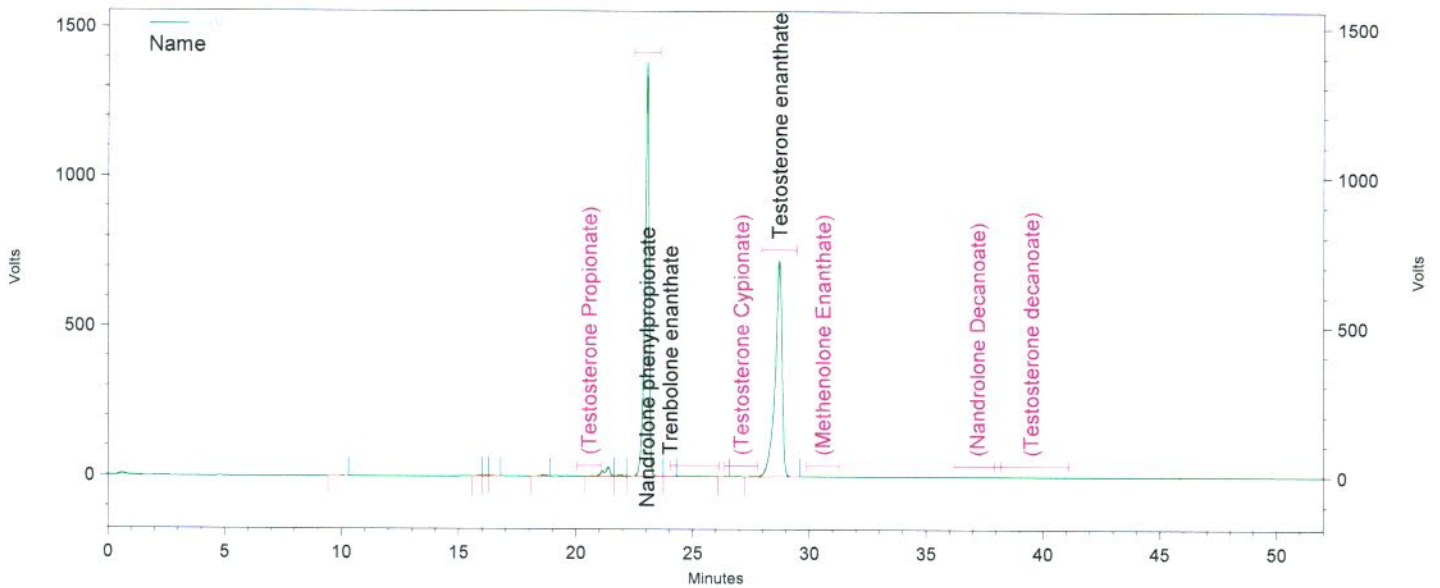

Channel A

| UV Results                  |        |          |           |              |
|-----------------------------|--------|----------|-----------|--------------|
| Name                        | Time   | Area     | Amount    | Area Percent |
|                             | 10.020 | 55202    | 0.000     | 0.042        |
|                             | 15.880 | 110852   | 0.000     | 0.084        |
|                             | 16.153 | 124429   | 0.000     | 0.095        |
|                             | 16.383 | 100980   | 0.000     | 0.077        |
|                             | 18.623 | 243060   | 0.000     | 0.185        |
| Testosterone Propionate     |        |          | 0.000 BDL |              |
|                             | 21.397 | 2199944  | 0.000     | 1.675        |
|                             | 21.917 | 250548   | 0.000     | 0.191        |
| Nandrolone phenylpropionate | 23.047 | 68727693 | 0.000     | 52.325       |
| Trenbolone enanthate        | 24.027 | 48360    | 0.000     | 0.037        |
|                             | 26.360 | 47327    | 0.000     | 0.036        |
| Testosterone Cypionate      |        |          | 0.000 BDL |              |
| Testosterone enanthate      | 28.710 | 59439365 | 0.000     | 45.253       |
| Methenolone Enanthate       |        |          | 0.000 BDL |              |


Nandrolone 0.000 BDL  
Decanoate  
Testosterone 0.000 BDL  
decanoate

|        |  |           |       |         |
|--------|--|-----------|-------|---------|
| Totals |  | 131347760 | 0.000 | 100.000 |
|--------|--|-----------|-------|---------|

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\November 2016\09-11-2016\Diverse  
Steroide\005\_STD Nandrolone Phenylprop und Test Enanthate-Rep5.dat  
Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Steroide A.met

Report steroid analysis

HPLC: Nr. 6  
 Analysendatum: 10.11.2016 05:03:26 Druckdatum: 14.11.2016 15:11:18  
 Vial: 3 Injektvolumen: 20  
 Probenidentifikation: STD Nandrolone Phenylprop und Test Enanthate  
 Operator: admin  
 Chrom.-Beding. ACN / Wasser, gradient  
 Säule Nr. 119, 100-5 C18 250 x 4.6 mm  
 Detektion: UV 240nm

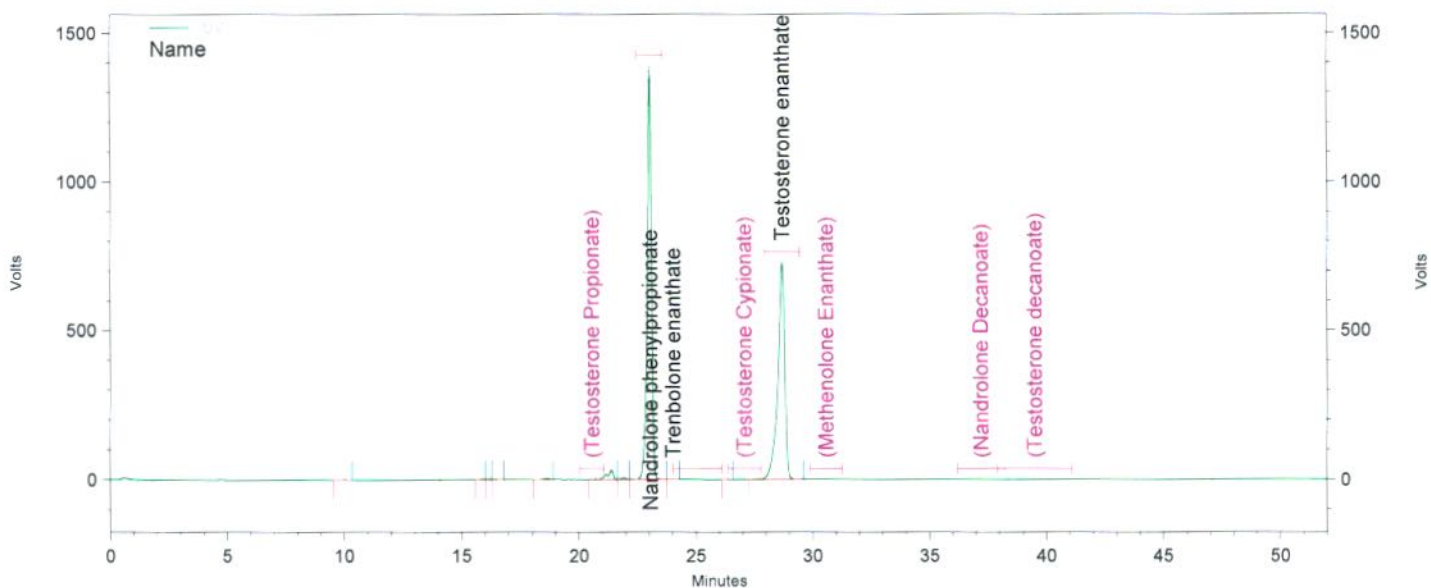

Channel A

| UV Results                  |        |          |           |              |
|-----------------------------|--------|----------|-----------|--------------|
| Name                        | Time   | Area     | Amount    | Area Percent |
|                             | 10.023 | 55853    | 0.000     | 0.042        |
|                             | 15.880 | 112453   | 0.000     | 0.085        |
|                             | 16.153 | 125271   | 0.000     | 0.095        |
|                             | 16.387 | 102659   | 0.000     | 0.078        |
|                             | 18.627 | 245098   | 0.000     | 0.186        |
| Testosterone Propionate     |        |          | 0.000 BDL |              |
|                             | 21.400 | 2206562  | 0.000     | 1.672        |
|                             | 21.923 | 251253   | 0.000     | 0.190        |
| Nandrolone phenylpropionate | 23.053 | 69064044 | 0.000     | 52.334       |
| Trenbolone enanthate        | 24.027 | 48080    | 0.000     | 0.036        |
|                             | 26.360 | 44960    | 0.000     | 0.034        |
| Testosterone Cypionate      |        |          | 0.000 BDL |              |
| Testosterone enanthate      | 28.710 | 59710756 | 0.000     | 45.247       |
| Methenolone Enanthate       |        |          | 0.000 BDL |              |

Nandrolone 0.000 BDL  
Decanoate  
Testosterone 0.000 BDL  
decanoate

|        |  |           |       |         |
|--------|--|-----------|-------|---------|
| Totals |  | 131966989 | 0.000 | 100.000 |
|--------|--|-----------|-------|---------|

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\November 2016\09-11-2016\Diverse  
Steroide\005\_STD Nandrolone Phenylprop und Test Enanthate-Rep6.dat  
Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Steroide A.met