

# Anabolic Lab™

Analytical Lab Testing for  
Anabolic Steroid Harm Reduction

<https://anaboliclab.com>

Lab Report – September 29, 2016

Primo of the Gods  
Nandrolone Phenylprop (100 mg/mL)  
Nandrolone Phenylpropionate



If you like the attached report, please visit [our website](https://anaboliclab.com) to find out how you can support independent and objective lab testing.

# Test report

Number **C 47303**

**Object of analysis:** Nandrolone Phenylprop 100 mg/mL (Primo of the Gods)

**Customer:** AnabolicLab.com **Samples:** 1

**Batch:** - **Receipt:** 15.08.2016

**Subject:** Assay, Identification Nandrolone Phenylpropionate

**Method(s):** HPLC, SOP M 006

Sample / Analysis:	Method	Status	Result	Specification	complies
<b><u>Nandrolone Phenylpropionate:</u></b>					
Identification	HPLC	N	complies	complies	i.O.
Content	HPLC	N	101.3 mg/mL	100 mg/mL	-
TAMC	SOP M 006	N	< 1 CFU/mL	-	-
TYMC	SOP M 006	N	< 1 CFU/mL	-	-

Method-status: G = GMP A = accredited V = generally validated P = validated on product N = not validated  
 E = external lab

**Remarks:** The sample(s) mentioned above have been analysed as they have been sent to us by the client. We have not controlled external processes like manufacturing, labelling, sampling, shipping, storage. These results are for information only and do not compensate for correct quality control by the manufacturer/distributor. Generally, pharmaceutical products have to be produced and distributed under full GMP/GDP regime, including GMP-compliant analyses with validated methods.

**Signatures:**  **controlled & released:**  **Staff:**


**completed:**  
29.09.2016

The enclosed results refer exclusively to the object of examination described above. The measuring accuracy is available on request. Only handwritten signatures are valid. Without written permission, it is not allowed to publish single parts of this report.

Report steroid analysis

HPLC: Nr.5  
 Analysendatum: 29.09.2016 00:29:51  
 Vial: 20  
 Probenidentifikation: C47303 A  
 Operator: admin  
 Chrom.-Beding.: ACN / Wasser, gradient  
 Temperatur: 40 °C  
 Säule: Nr. 135, 250 x 4,6, C18  
 Detektion: UV 240nm

Druckdatum: 29.09.2016 09:31:55  
 Injektvolumen: 20

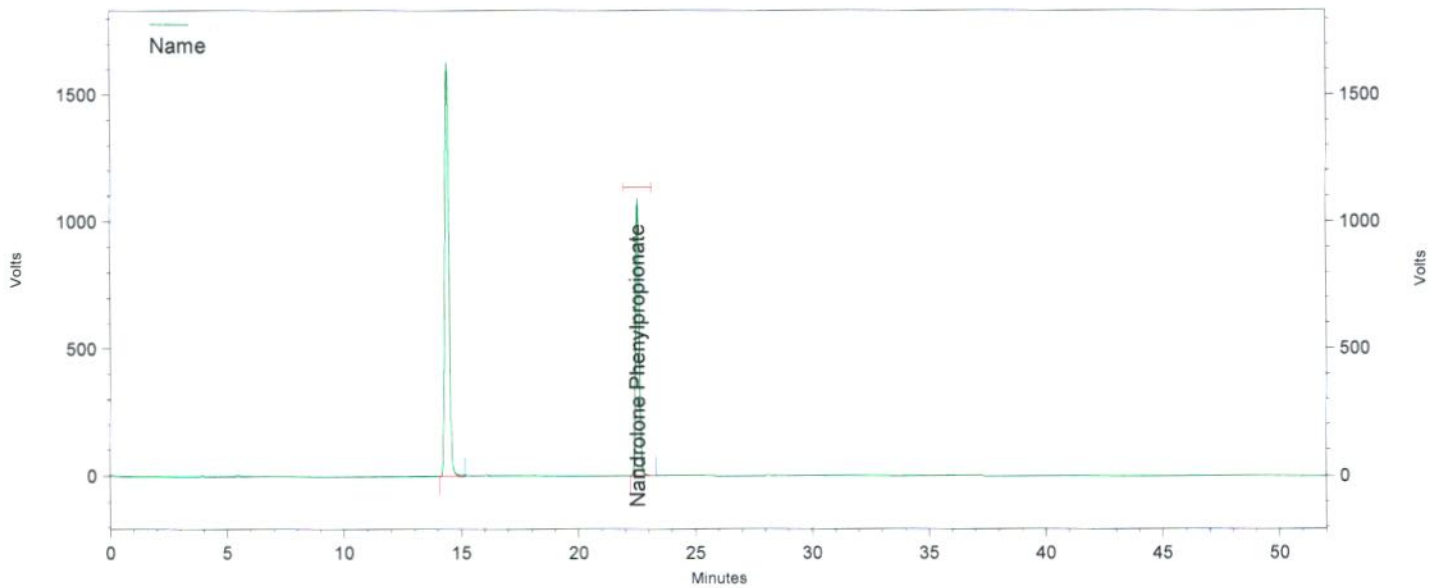

Channel A


UV Results				
Name	Time	Area	Amount	Area Percent
	14,397	78962608	0,000	63,637
Nandrolone Phenylpropionate	22,530	45119518	0,000	36,363
Totals		124082126	0,000	100,000

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\September 2016\16-09-26  
 Steroide\026\_C47303 A-Rep1.dat  
 Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Steroide A.met

Report steroid analysis

HPLC: Nr.5  
 Analysendatum: 29.09.2016 01:23:20  
 Vial: 20  
 Probenidentifikation: C47303 A  
 Operator: admin  
 Chrom.-Beding.: ACN / Wasser, gradient  
 Temperatur: 40 °C  
 Säule: Nr. 135, 250 x 4,6, C18  
 Detektion: UV 240nm

Druckdatum: 29.09.2016 09:32:00  
 Injektvolumen: 20

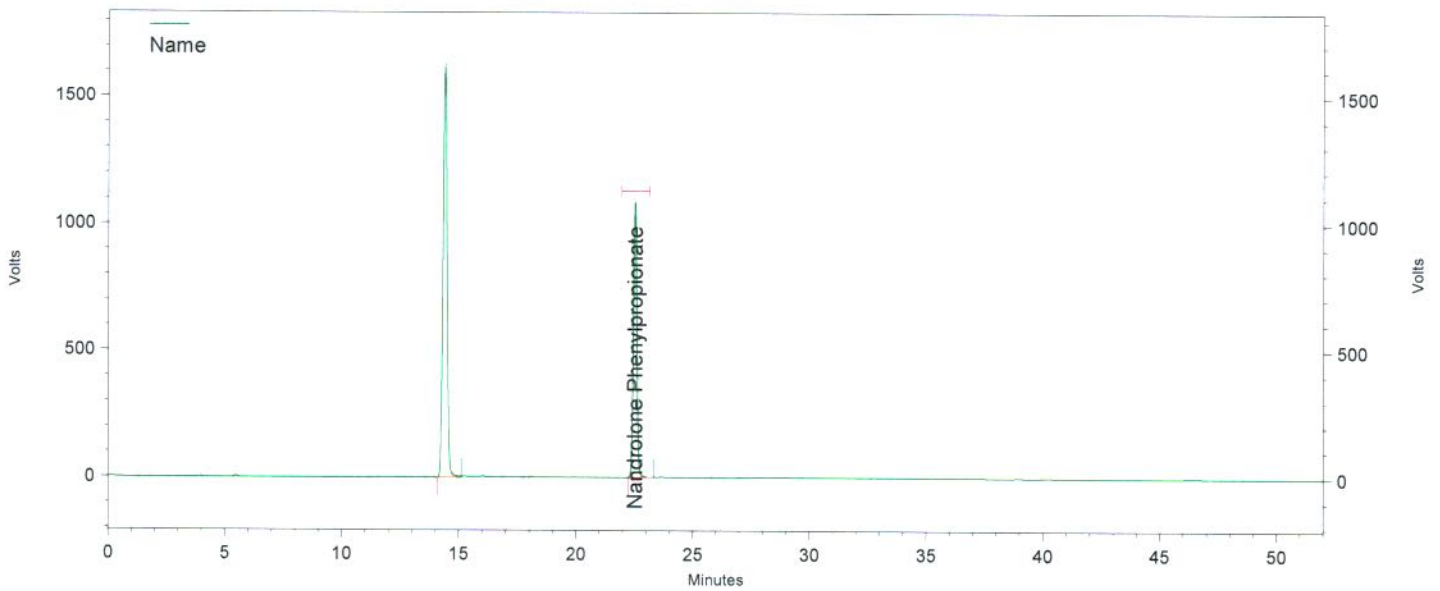

Channel A

UV Results

Name	Time	Area	Amount	Area Percent
Nandrolone Phenylpropionate	14,390	79039298	0,000	63,655
	22,523	45128213	0,000	36,345


Totals		124167511	0,000	100,000
--------	--	-----------	-------	---------

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\September 2016\16-09-26  
 Steroide\026\_C47303 A-Rep2.dat  
 Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Steroide A.met

Report steroid analysis

HPLC: Nr.5  
 Analysendatum: 29.09.2016 02:16:45  
 Vial: 20  
 Probenidentifikation: C47303 A  
 Operator: admin  
 Chrom.-Beding. ACN / Wasser, gradient  
 Temperatur: 40 °C  
 Säule Nr. 135, 250 x 4,6, C18  
 Detektion: UV 240nm

Druckdatum: 29.09.2016 09:32:04  
 Injektvolumen: 20

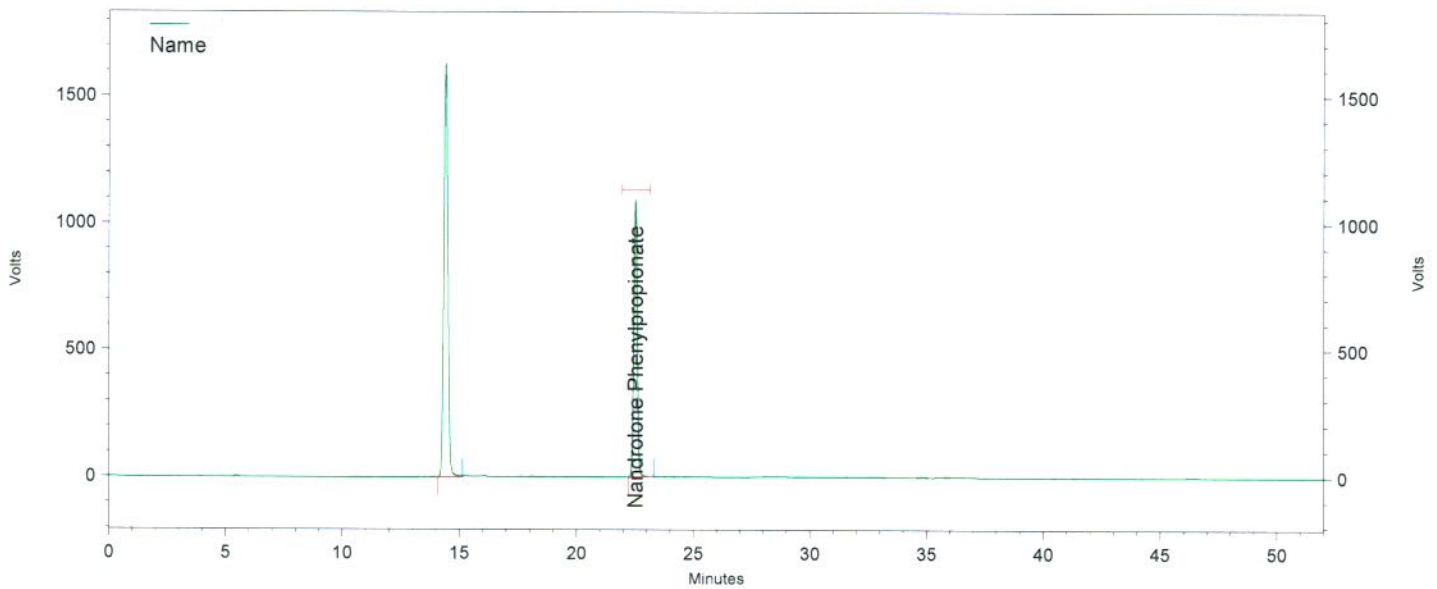

Channel A

UV Results

Name	Time	Area	Amount	Area Percent
Nandrolone Phenylpropionate	14,407	79258441	0,000	63,634
	22,520	45294806	0,000	36,366


Totals		124553247	0,000	100,000
--------	--	-----------	-------	---------

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\September 2016\16-09-26  
 Steroide\026\_C47303 A-Rep3.dat  
 Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Steroide A.met

Report steroid analysis

HPLC: Nr.5  
 Analysendatum: 29.09.2016 03:10:10  
 Vial: 21  
 Probenidentifikation: C47303 B  
 Operator: admin  
 Chrom.-Beding. ACN / Wasser, gradient  
 Temperatur: 40 °C  
 Säule Nr. 135, 250 x 4,6, C18  
 Detektion: UV 240nm

Druckdatum: 29.09.2016 09:32:08  
 Injektvolumen: 20

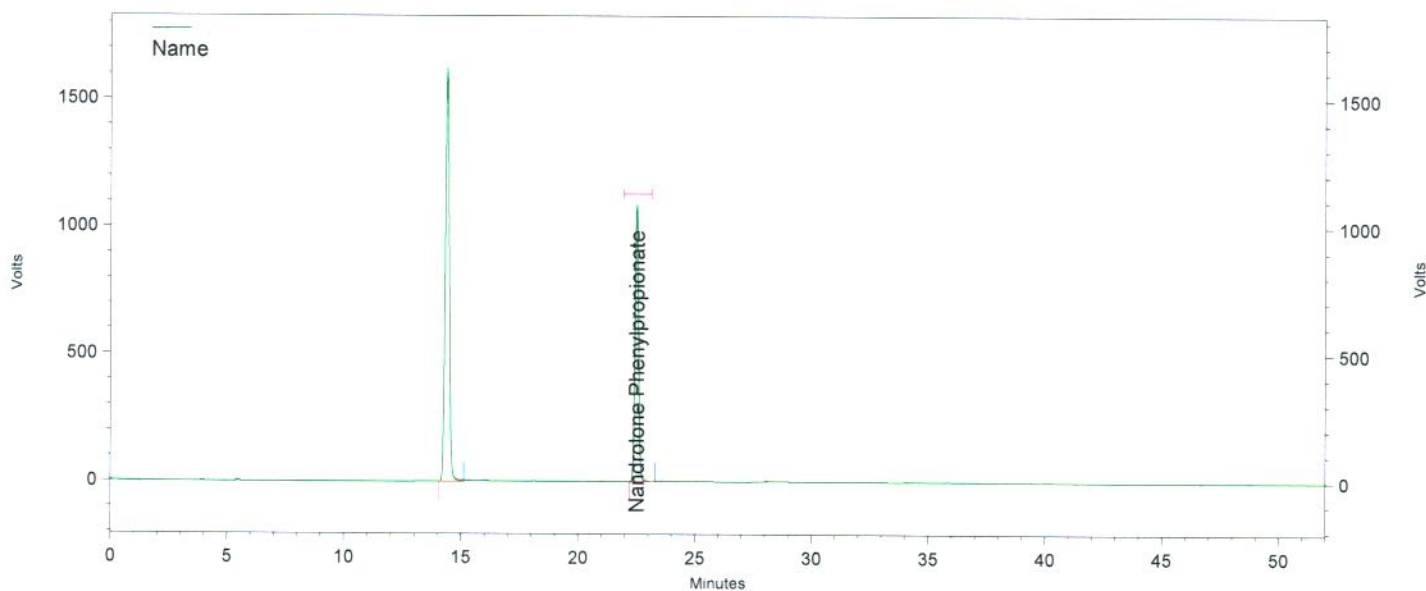

Channel A

UV Results

Name	Time	Area	Amount	Area Percent
Nandrolone Phenylpropionate	14,393	79002302	0,000	63,667
	22,507	45085334	0,000	36,333


Totals		124087636	0,000	100,000
--------	--	-----------	-------	---------

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\September 2016\16-09-26  
 Steroide\027\_C47303 B-Rep1.dat  
 Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Steroide A.met

Report steroid analysis

HPLC: Nr.5  
 Analysendatum: 29.09.2016 04:03:36  
 Vial: 21  
 Probenidentifikation: C47303 B  
 Operator: admin  
 Chrom.-Beding. ACN / Wasser, gradient  
 Temperatur: 40 °C  
 Säule Nr. 135, 250 x 4,6, C18  
 Detektion: UV 240nm

Druckdatum: 29.09.2016 09:32:13  
 Injektvolumen: 20

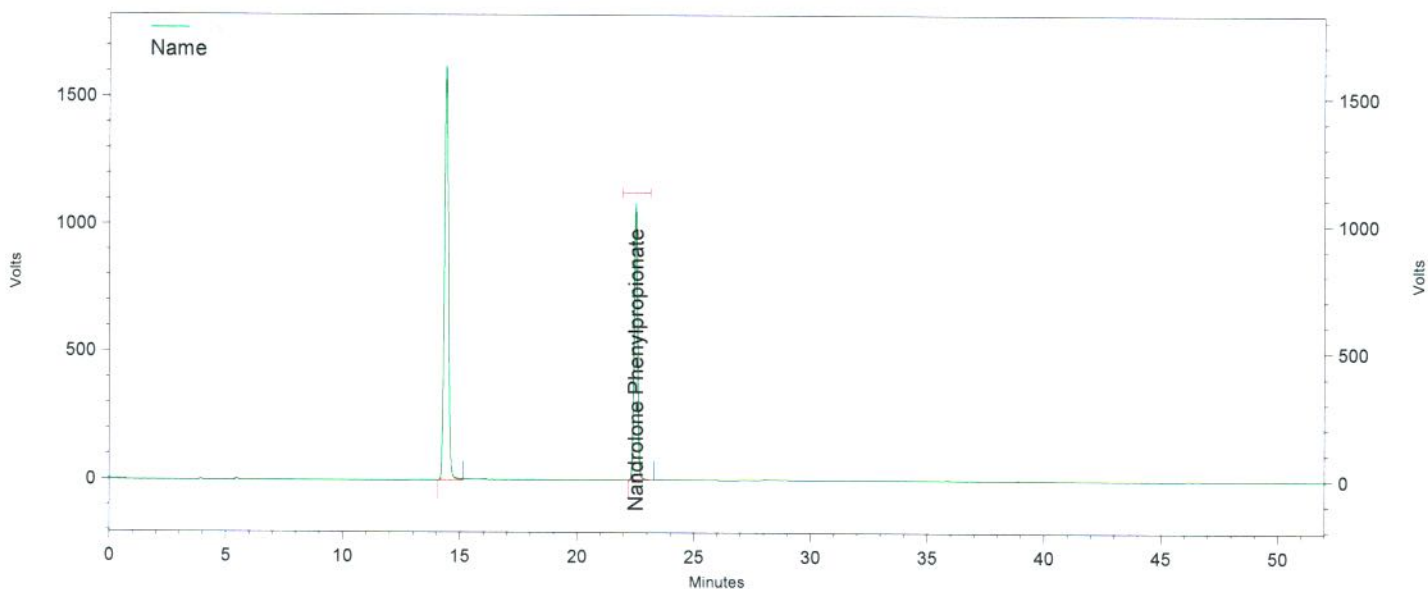

Channel A

UV Results

Name	Time	Area	Amount	Area Percent
Nandrolone Phenylpropionate	14,403	78963910	0,000	63,690
	22,513	45018001	0,000	36,310


Totals		123981911	0,000	100,000
--------	--	-----------	-------	---------

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\September 2016\16-09-26  
 Steroide\027\_C47303 B-Rep2.dat  
 Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Steroide A.met

Report steroid analysis

HPLC: Nr.5  
 Analysendatum: 29.09.2016 04:57:00  
 Vial: 21  
 Probenidentifikation: C47303 B  
 Operator: admin  
 Chrom.-Beding. ACN / Wasser, gradient  
 Temperatur: 40 °C  
 Säule Nr. 135, 250 x 4,6, C18  
 Detektion: UV 240nm

Druckdatum: 29.09.2016 09:32:17  
 Injektvolumen: 20

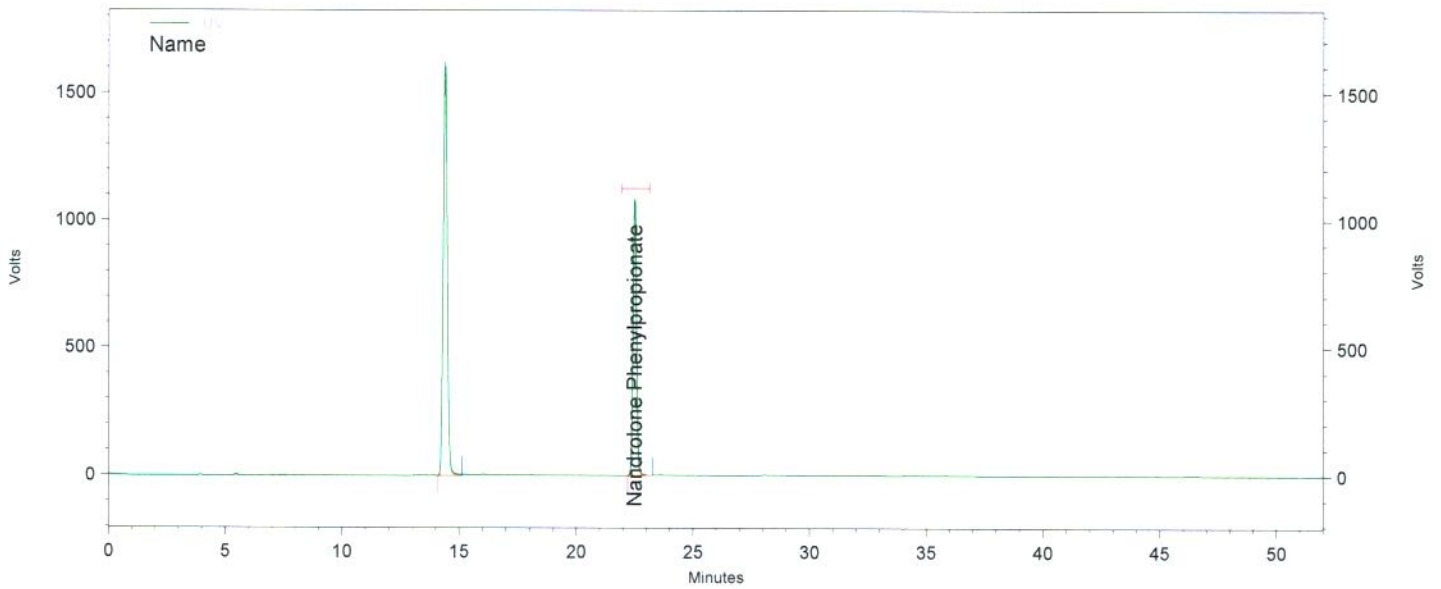

Channel A

UV Results

Name	Time	Area	Amount	Area Percent
Nandrolone Phenylpropionate	14,393	79144732	0,000	63,696
	22,503	45109103	0,000	36,304

Totals		124253835	0,000	100,000
--------	--	-----------	-------	---------

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\September 2016\16-09-26  
 Steroide\027\_C47303 B-Rep3.dat  
 Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Steroide A.met